


# Paraffin Block Mailer

## M477-6

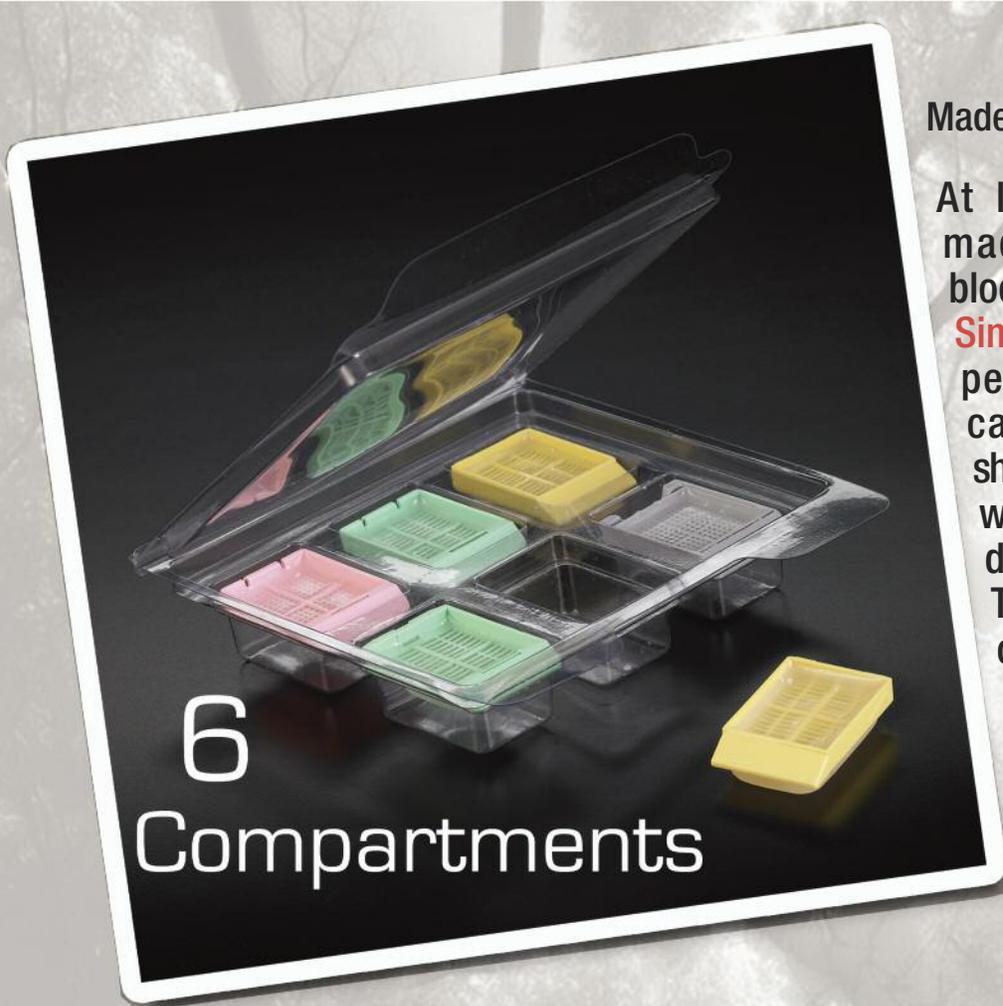

Made of PVC

At last, a transport container made especially for paraffin blocks. Sturdy and easy to use, the **Simport** mailer is transparent for perfect viewing of contents. It can be used for handling and shipping up to 6 blocks as needed, while having to keep only one model of block mailer in inventory. The attached cover is easy to close but very secure when closed. Suitable for all regular models of tissue and biopsy cassettes.

**Cat. #**

M477-6

**Dimensions**

134 x 141 x 29 mm H  
(5 1/4 x 5 5/8 x 1 1/8 in. H)

**Qty/Cs**



50

# Simport

[www.simport.com](http://www.simport.com)



# THE DispoCut™ DISSECTING BOARD M630

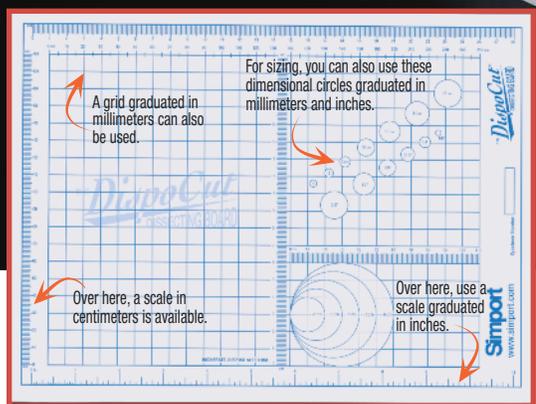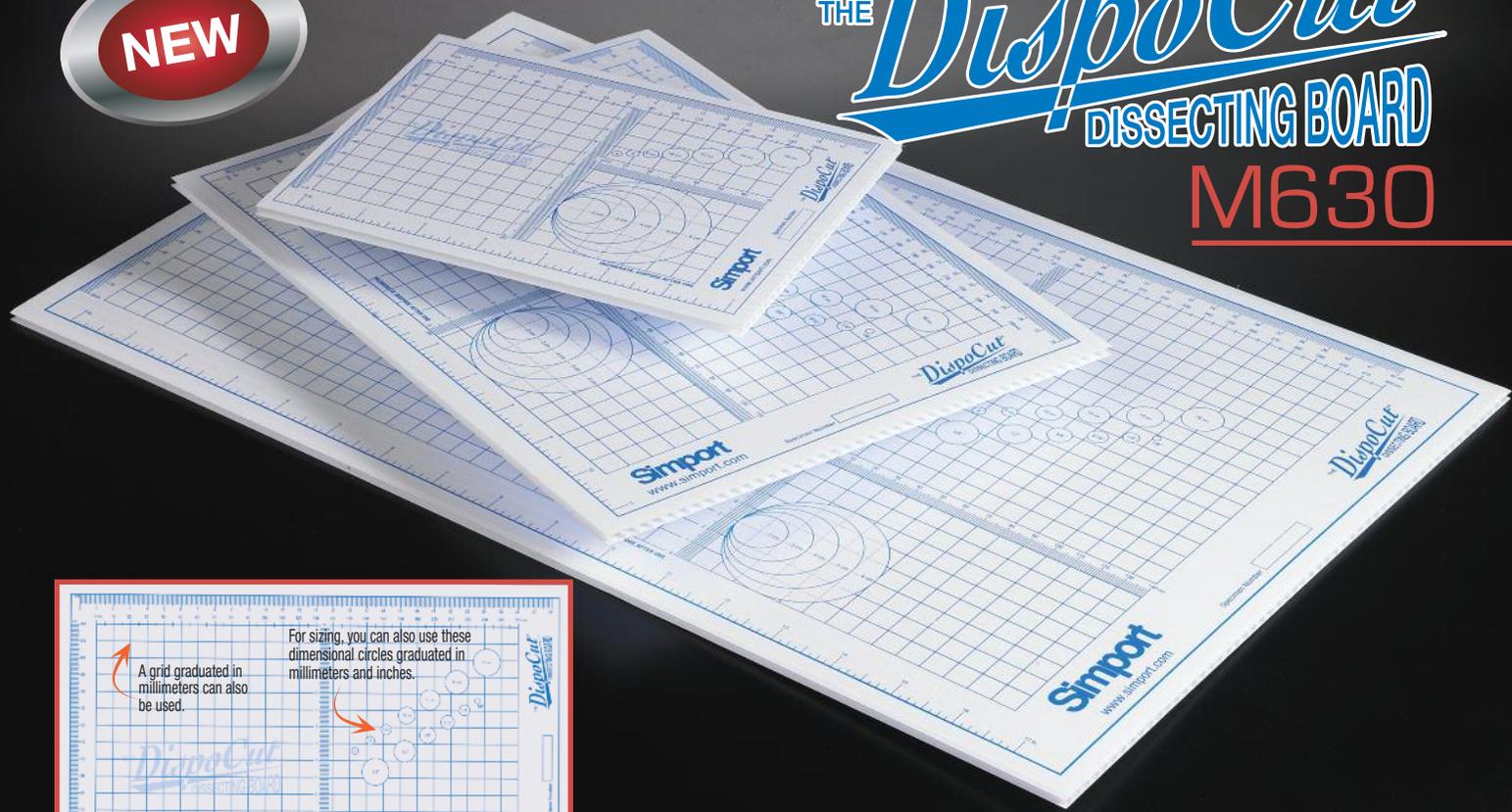

| Cat. # | Dimensions                 | Qty/Pk | Qty/Cs |
|--------|----------------------------|--------|--------|
| M630-1 | 6 x 8 in. (152 x 203 mm)   | 24     | 96     |
| M630-2 | 9 x 12 in. (229 x 305 mm)  | 12     | 48     |
| M630-3 | 12 x 19 in. (305 x 483 mm) | 12     | 24     |

## M630 Disposable Dissecting Board

The DispoCut™ Disposable Dissecting Board is strong yet inexpensive enough to throw away. It is especially developed to provide a clean, safer and more efficient way of handling infectious tissue specimens for the pathologist. It can be used on both sides, a great money saving feature. Printed with helpful imperial and metric dimensional scales in inches, centimeters and millimeters. Available in three sizes to accommodate small biopsies to large gross anatomy procedures.

DispoCut™ is extrusion made using a copolymer resin in order to increase impact and high and low temperature performance. Copolymer resins are also used because they retain the ability to be flexed an unlimited number of times without breaking. Chemically, DispoCut is inert. At regular temperatures most oils, solvents and water have no effect, allowing it to perform under adverse weather conditions or as a product component exposed to harsh chemicals.

Each surface has many conveniently located rulers including seven dimensional circles from 4 to 16 mm in diameter and five larger ones from 2 to 6 cm in diameter. Surface will not dull knives. Perfect for use with dissecting pins. All sizes of DispoCut™ Disposable Dissecting Boards are ideal for pathology labs, medical schools, classroom dissections, forensic labs, medical research labs, pharmaceutical labs, surgery and more.

Temperature performance range: -17 °F to 160 °F, Melting point: 162 °C, 324 °F

Sterilization: DispoCut™ may be wiped down with alcohol with no effect on the material. There are many types of autoclaves on the market. DispoCut™ is being autoclaved (in excess of 3 years) in many situations. There are some extreme settings which will be damaging and we suggest a trial to determine performance in your equipment. Ethylene oxide gas sterilization has no effect of any kind on DispoCut™ and may be repeated any number of times.

Physical Data: Specific Gravity 0.90 - 0.96. Essentially insoluble in water

For further details, contact our Customer Service Department at: (450) 464-1723

DispoCut is a trademark of Simport Scientific.



2588 Bernard-Pilon  
Beloil, QC J3G 4S5 Canada

Tel: (450) 464-1723  
Fax: (450) 464-3394

E-mail: info@simport.com  
Website: www.simport.com

Distributed by: