

Simport[®] LABORATORY

Canadian and Family Owned Since 1975

EcoTainer 24[™]



Degradable is now available for the laboratory. Make it your choice and protect your environment.



The Simport EcoTainer 24[™] will biodegrade to become some of the soil's organic components in less than 7 years instead of up to 400 years when using conventional plastics.

www.simport.com



INTRODUCING THE NEXT GENERATION OF 24-HR URINE COLLECTION CONTAINERS



The Simport® URISAFE®, one of the most innovative 24-HR Urine Container on the market today, is now available in a new version: The EcoTainer 24™, a rapidly degradable URISAFE® Container.

Plastics take hundreds of years to degrade naturally in the environment. The Simport® EcoTainer 24™, will biodegrade to become some of the soil's organic components in less than 7 years instead of up to 400 years when using conventional plastics. The only condition that is necessary for the EcoTainer 24™ to biodegrade is constant contact with other degrading material. No heat, physical stress, oxygen or sunlight necessary.

The EcoTainer 24™ has the ability to break down safely and quickly, by biological means, into the raw materials of nature and disappear into the environment.


For IVD use

- Does not contain heavy metals
- No negative effect to its physical properties
- No heat, light or mechanical stress needed for product breakdown
- No toxic residue
- Fully degradable within 7 years, aerobically or anaerobically

Image of high-density polyethylene surface after 12 months soil burial



Non biodegradable surface

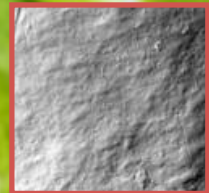

Surface of the EcoTainer 24™

Uniquely designed to be the most user friendly



- 1 The unique snap valve pour spout is easily popped open and offers dripless pouring
- 2 Specially designed leakproof screw cap with liner for safer transport
- 3 Large centrally located handle for good balance and ease of manipulation
- 4 In regular position, sample volume can be read in 100 ml increments
- 5 Place container in upright position and read in 50 ml increments

The unique snap valve pour spout:

- Controlled flow rate for better handling of poured volumes
- Reduced risks of aerosol contamination when pouring
- Eliminates splashing and exposure to hazardous body fluids

Low form design:

- Convenient refrigerator storage with less wasted space
- Anatomically designed for ease of patient use
- Large 69 mm screw lid

Chemical and physical resistance

- Can be subjected to freezing, thawing and EtO gas sterilization without causing changes in materials or physical appearance.
- Metal free and resistant to hydrochloric acid

Cat. #	Vol. (Liters)	Dimensions (cm)	Qty/Cs
B350ECO	3	11.5 x 24.5 x 16.0 H	40

EcoTainer 24 is a trademark of Simport® Scientific.
Urisafe is a registered trademark of Simport® Scientific.

Simport[®]
Scientific

2588 Bernard-Pilon
Beloil, Qc J3G 4S5 Canada

Phone: (450) 464-1723

Fax: (450) 464-3394

Email: info@simport.com

Website: www.simport.com

Distributed by: