


# Simport®

*Canadian and Family Owned Since 1975*

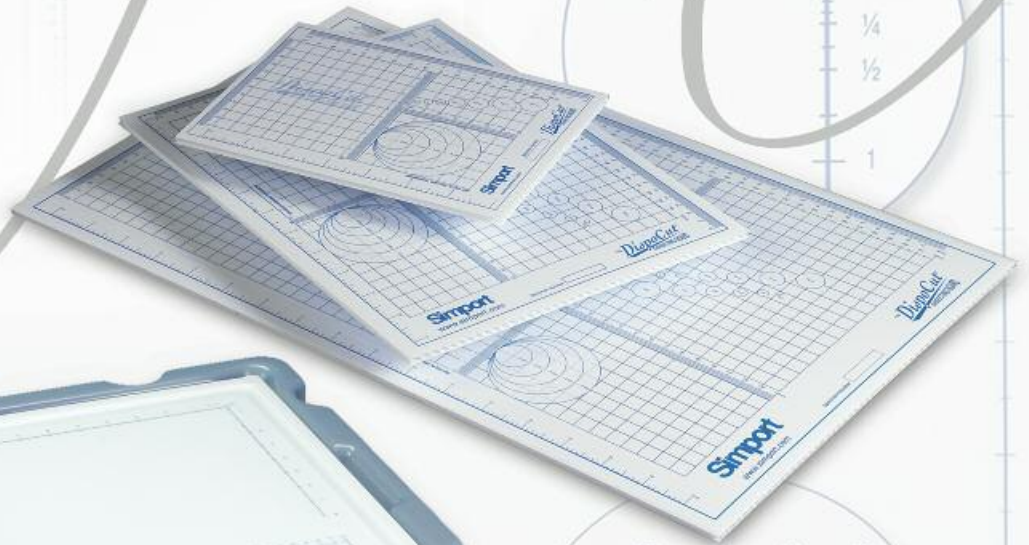

**The First and Only  
Reversible Dissecting Table  
for the Laboratory**


[www.simport.com](http://www.simport.com)




# The DissectTable™

## Simpport


Flip it and have access to two smaller tables also supplied with dimensional scales.

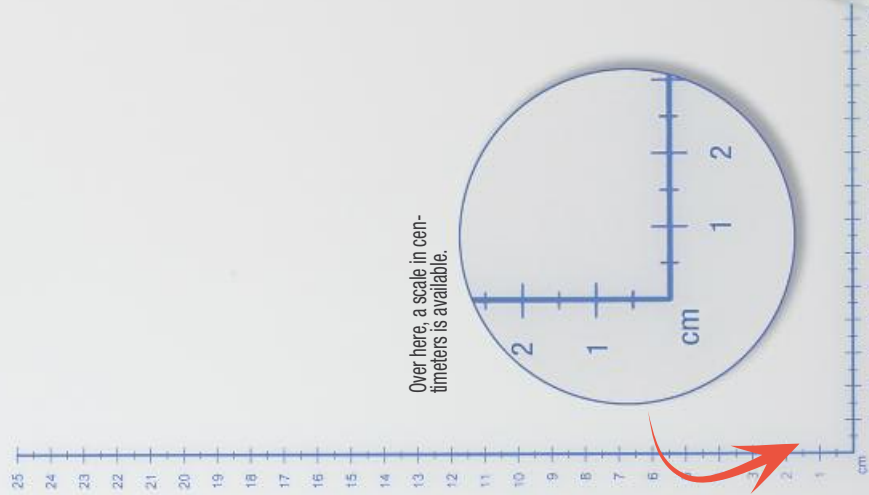


Over here, a scale in centimeters is available.

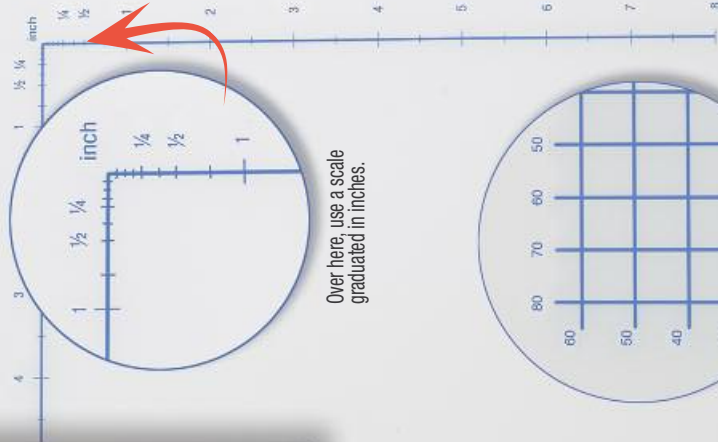


Over here, use a scale graduated in inches.

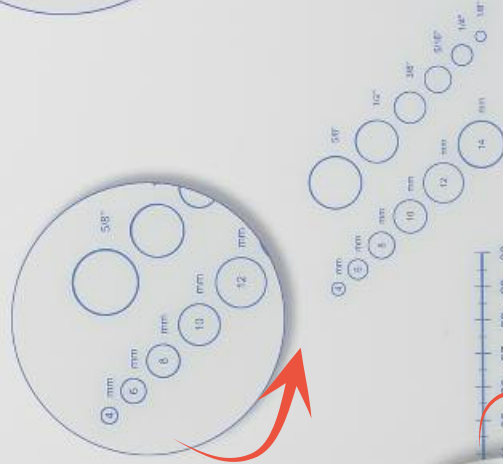

For sizing, you can also use these dimensional circles graduated in millimeters and inches.

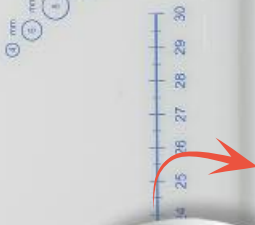

A drain groove is carved all around the edge.

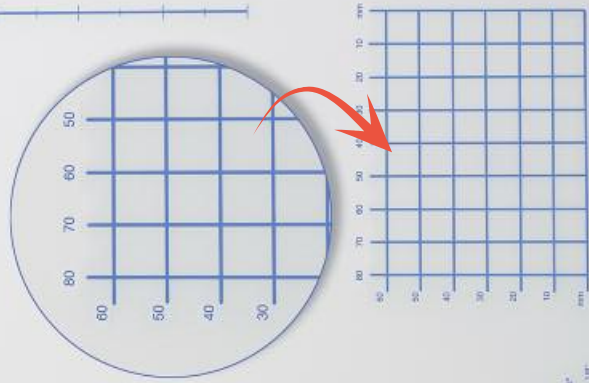

A grid graduated in millimeters can also be used.



# M620


## DissecTable™ Dissecting Board

Made of high-density polyethylene

- Large working area on one side, two smaller ones on the other side
- Will not dull fine surgical blades
- Carved drain groove all around the edge
- Dimensional scales in centimetres and inches
- 60 x 80 mm grid made of forty-eight 10 mm squares
- Six dimensional circles with a diameter from 4 to 14 mm and 1/8 to 5/8 in.
- Rubber feet under corners ensuring more stability

A new and unique approach makes this dissecting board more convenient than any other found on the market today. Made of heavy-duty stain resistant thick polyethylene, it will last for years to come without, bending, swelling or changing shape.

Dimensions: 575 mm x 400 mm x 12.5 mm (23 x 16 x 1/2 in.)

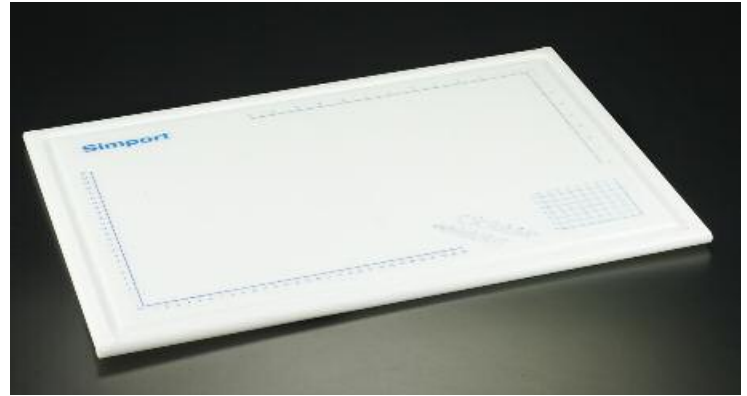

Cat. #	Description	Qty/Cs
M620	DissecTable™	1

# M618

## DissecTable™ Jr. Dissecting Board

Made of high-density polyethylene

A smaller DissecTable™ is also available with the same features and benefits as the M620. Perfect for smaller counter areas. Dimensions: 330.2 mm x 279.4 mm x 12.5 mm (13 x 11 x 1/2 in H)

Cat. #	Description	Qty/Cs
M618	DissecTable™ Jr.	1




# M625 DissecTable™ Board Base

Made of high-density polyethylene

To make dissecting more comfortable, this heavy-duty base is used to elevate the M620 DissecTable™ Board to the right height. The bases are stackable and will not move sideways during the dissecting work. The base will also retain excess fluid if necessary.

Dimensions: 481 mm x 656 mm x 91 mm (19 1/4 x 26 1/4 x 3 5/8 in H)

Cat. #	Description	Qty/Cs
M625	DissecTable™ Board Base	1

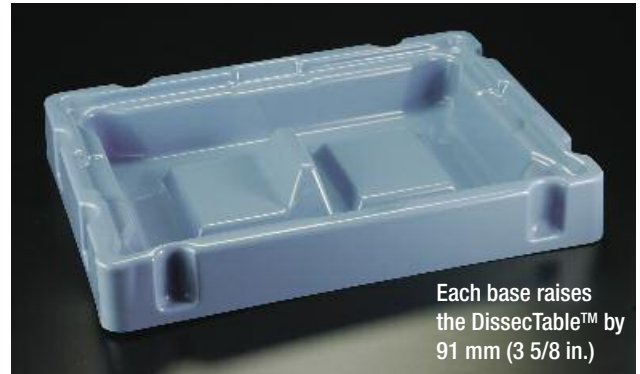


For better comfort, use stable base to hold The DissecTable™ Board at the proper height.



To elevate the DissecTable™ Board even more, you can add as many board bases as needed.

# THE *DispoCut*<sup>TM</sup> DISSECTING BOARD

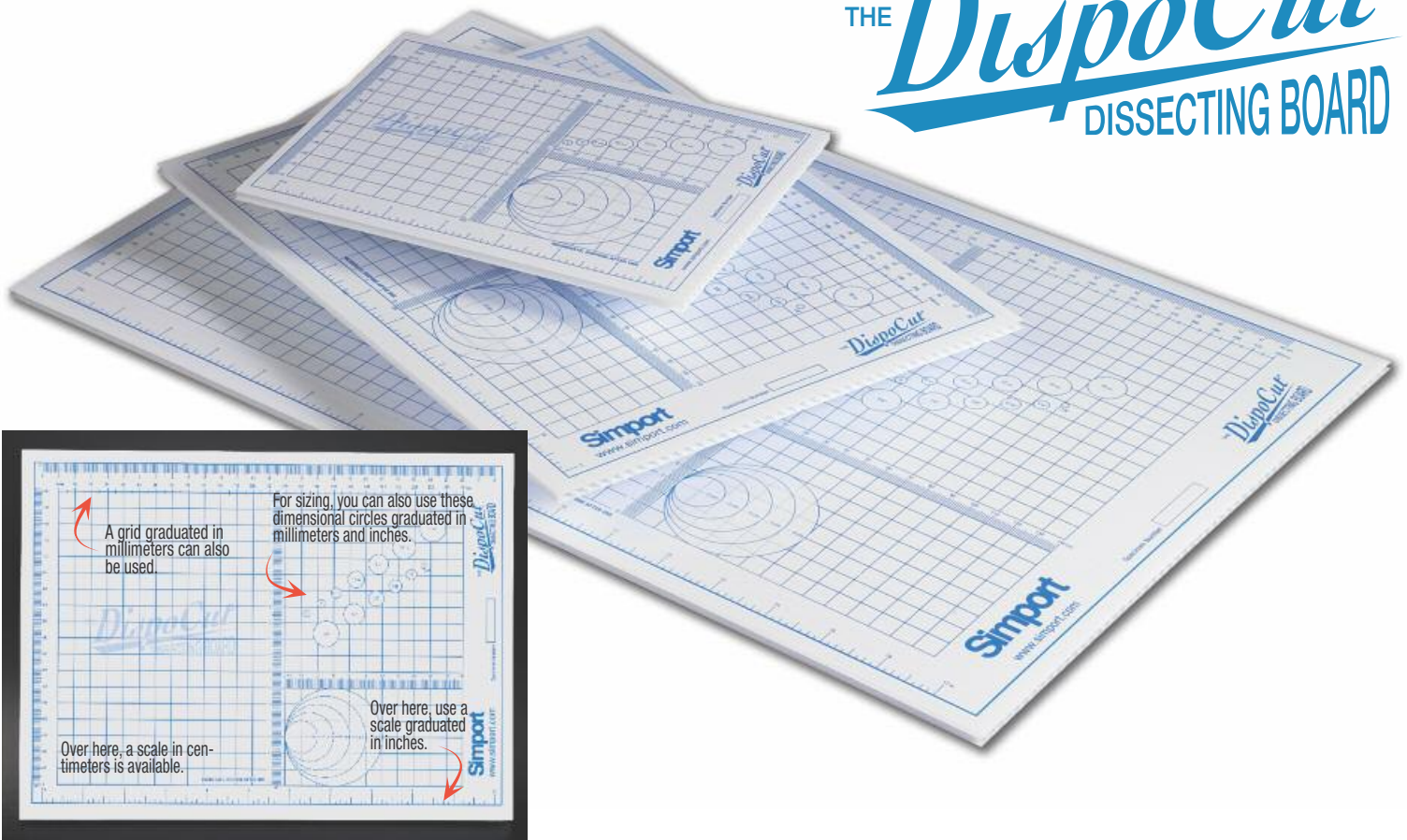

## M630 Disposable Dissecting Board

- Reversible
- Printed imperial and metric dimensional scales in inches, centimeters and millimeters
- Seven dimensional circles from 4 to 16 mm and five larger ones from 2 to 6 cm in diameter
- Strong surface will not dull knives but well suited for using dissecting pins
- Temperature performance range: -27 °C to 71 °C (17 °F to 160 °F)
- Melting point: 162 °C (324 °F)

The DispoCut<sup>TM</sup> Disposable Dissecting Board is strong yet inexpensive enough to throw away. It is especially developed to provide a clean, safer and more efficient way of handling infectious tissue specimens. Ideal for pathology labs, medical schools, classroom dissections, forensic labs, medical research labs, pharmaceutical labs, surgery and more. When sterilizing, DispoCut<sup>TM</sup> may be safely wiped down with alcohol. DispoCut<sup>TM</sup> is being autoclaved in many situations. Some other cleaning agents might be damaging and we suggest a trial to determine performance. Ethylene oxide gas sterilization has no negative effect on DispoCut<sup>TM</sup> and may be repeated any number of times.

Available in three sizes.

Cat. #	Dimensions	Qty/Pk	Qty/Cs
M630-1	152 x 203 mm (6 x 8 in.)	24	96
M630-2	229 x 305 mm (9 x 12 in.)	12	48
M630-3	305 x 483 mm (12 x 19 in.)	12	24

DissecTable and Dispocut are trademark of Simport® Scientific

**Simport**<sup>®</sup>  
*Scientific*

2588 Bernard-Pilon  
Beloeil, Qc J3G 4S5 Canada

Phone: (450) 464-1723

Fax: (450) 464-3394

Email: [info@simport.com](mailto:info@simport.com)

Website: [www.simport.com](http://www.simport.com)

Distributed by: