

Simport®

Canadian and Family Owned Since 1975

CoreDish™

Multiple Biopsy Containers

Half prefilled with
10% formalin

Also available
without formalin

If you **REALLY** care about
your biopsies, let us help you
FACILITATE their handling!



www.simport.com

CoreDish™

MULTIPLE BIOPSY CONTAINERS HALF PREFILLED WITH 10% FORMALIN

Made of polystyrene

Few recommendations concerning how the biopsies should be handled have been published. Performing a large number of biopsies means an increase in the number of containers handled and consequently a technical overload of the transmission network, which occurs without any financial counterpart. A new approach had to be developed in order to increase productivity.

Simport® is proud to offer a multi-compartment container in the shape of a dish and half prefilled with 10% Neutral Buffered Formalin, for holding and transporting biopsies. It is supplied with a leakproof closure with O-ring ensuring total protection of contents. It conforms to OSHA directives. **The Simport® CoreDish™ measures only 15 x 95 mm in diameter.** Each compartment is clearly identified to allow proper placement and visualization of the biopsy being inserted. **Thanks to the CoreDish™ it is no more necessary to use a multitude of individual containers, thereby reducing risks of confusion.**

The Simport® CoreDish™ offers many configurations in order to hold different biopsies of the breast, prostate, upper GI tract and lower GI tract. A label allows essential information to be written such as patient I.D., doctor, date and time.


The CoreDish™ is also available without formalin. See Series M971.

For IVD use 




The CoreDish uses up to 75% less formalin.


1. Well labelled and easy to read identification of biopsy origin
2. Each biopsy is isolated from each other by a divider
3. Molded ridges around lid make it easy to open and close

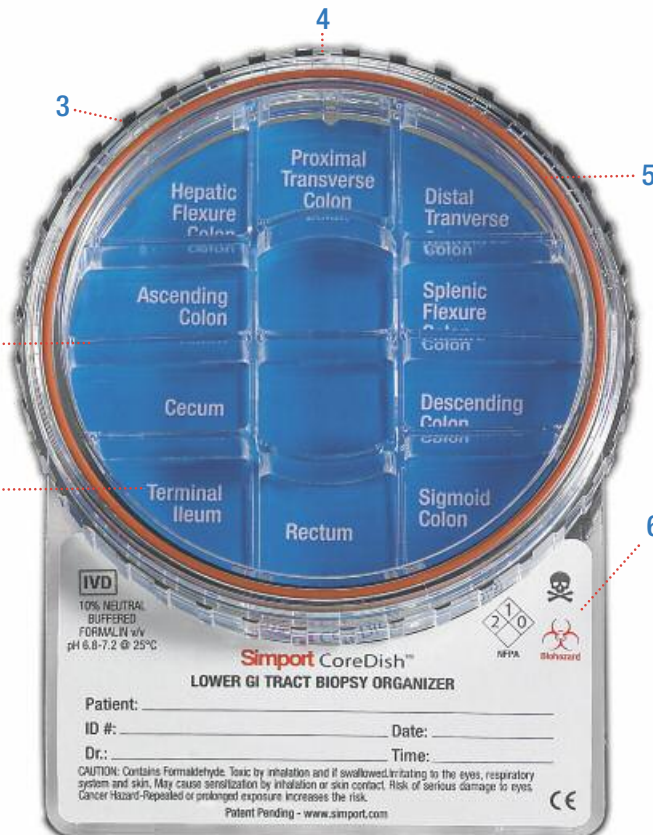

4. Fast thread design allows lid to be removed with only 1 1/2 turn
5. Specially designed O-ring for a more secure and positive leakproof seal
6. Warning label has space for patient identification

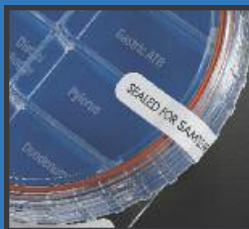

A security label is supplied to ensure integrity of contents from collection stage to re-opening of CoreDish™.

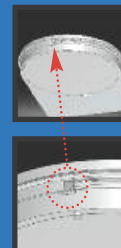

To close, turn clockwise until you feel a firm stop.



M970-D5B-1

BREAST BIOPSY CONTAINER For separation, imaging and transport of core needle specimens

Designed specifically for radiography, separation, imaging and transport of core needle breast biopsies. Special absorbent liner keeps specimens moist (when saline solution is added) prior to radiography while helping to attenuate the x-ray beam. Four compartments are clearly identified (3, 6, 9 and 12) and the radiolucent numbers show up clearly on the radiograph. Formalin may be added prior to transportation to pathology for analysis. An area for writing patient information is provided on the label. Leakproof seal, thanks to O-ring lid, allows for safe and easy transport of the specimens from collection to analysis. **Not available in USA.**

Cat. #	For	Compartments	Qty/Pk	Qty/Cs
M970-D5B-1	Breast	5	1	10



M970-D5B-2

BREAST BIOPSY CONTAINER Half Prefilled with 10% Formalin

Simport® is proud to offer a multi-compartment container (out of five compartments, four are labelled: Left Upper Quadrant, Right Upper Quadrant, Left Lower Quadrant, Right Lower Quadrant) in the shape of a dish and half prefilled with 10% Neutral Buffered Formalin, for holding and transporting biopsies. It conforms to OSHA directives. Each compartment is clearly identified to allow proper placement and visualization of the breast biopsy being inserted. A writing area for patient information is provided. Leakproof seal, thanks to O-ring lid, allows for safe and easy transport of the specimens from collection to analysis.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D5B-2	Breast	5	YES	1	10
M971-D5B-2	Breast	5	NO	1	10



M970-D8P

PROSTATE BIOPSY CONTAINER Half Prefilled with 10% Formalin

For prostate biopsies. Eight compartments. Leakproof seal thanks to O-ring lid. An area for patient information is provided. Six labeled compartments: Base, Lateral Base, Medium, Lateral Medium, Apex, Lateral Apex. Leakproof seal, thanks to O-ring lid, allows for safe and easy transport of the specimens from collection to analysis.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D8P	Prostate	8	YES	1	10
M971-D8P	Prostate	8	NO	1	10



M970-D8UGI

UPPER GI TRACK BIOPSY CONTAINER Half Prefilled with 10% Formalin

For upper GI track biopsies. Eight compartments. Leakproof seal thanks to O-ring lid. An area for patient information is provided. Seven labeled compartments: Gastric Card, Gastric Body, GE Junction, Gastric ATR, Distal Esophage, Pylorus, Duodenum. Leakproof seal, thanks to O-ring lid, allows for safe and easy transport of the specimens from collection to analysis.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D8UGI	Upper GI Tract	8	YES	1	10
M971-D8UGI	Upper GI Tract	8	NO	1	10



M970-D12LGI

LOWER GI TRACK BIOPSY CONTAINER Half Prefilled with 10% Formalin

For lower GI track biopsies. Twelve compartments. Leakproof seal thanks to O-ring lid. An area for patient information is provided. Ten labeled compartments: Proximal Flexure Colon, Hepatic Flexure Colon, Distal Transverse Colon, Ascending Colon, Splenic flexure Colon, Cecum, Descending Colon, Terminal Ileum, Rectum, Sigmoid Colon. Leakproof seal, thanks to O-ring lid, allows for safe and easy transport of the specimens from collection to analysis.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D12LGI	Lower GI Tract	12	YES	1	10
M971-D12LGI	Lower GI Tract	12	NO	1	10



M970-D12P

PROSTATE BIOPSY CONTAINER Half Prefilled with 10% Formalin

For prostate biopsies. Twelve compartments. Leakproof seal thanks to O-ring lid. An area for patient information is provided. Twelve labeled compartments: L Base, R Base, L Lateral Base, R Lateral Base, L Lateral Medium, L Medium, R Medium, R Lateral Medium, L Lateral Apex, L Apex, R Apex, R Lateral Apex. Leakproof seal, thanks to O-ring lid, allows for safe and easy transport of the specimens from collection to analysis.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D12P	Prostate	12	YES	1	10
M971-D12P	Prostate	12	NO	1	10

CoreDish™

GENERAL PURPOSE MULTIPLE BIOPSY CONTAINERS

95kPa tested



M970-D5

5-COMPARTMENT BIOPSY CONTAINER Half Prefilled with 10% Formalin

Made of polystyrene

Performing a large number of biopsies means an increase in the number of containers handled and consequently a technical overload of the transmission network, which occurs without any financial counterpart. A new approach had to be developed in order to increase productivity. Simport® is proud to offer a five-compartment container in the shape of a dish and half prefilled with 10% Neutral Buffered Formalin, for holding and transporting biopsies. It is supplied with a leakproof screw closure with O-ring ensuring total protection of contents. It conforms to OSHA directives. A writing area for patient information is provided. Compartments are numbered from 1 to 5.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D5	General purpose	5	YES	1	10
M971-D5	General purpose	5	NO	1	10



M970-D8

8-COMPARTMENT BIOPSY CONTAINER Half Prefilled with 10% Formalin

Made of polystyrene

Simport® is proud to offer a multi-compartment container for up to 8 biopsies. The screw on lid incorporates an O-ring in order to make it leakproof and protect its contents. The CoreDish™ conforms to OSHA directives. An area for patient information is provided on the label. Compartments are numbered from 1 to 8.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D8	General purpose	8	YES	1	10
M971-D8	General purpose	8	NO	1	10




M970-D12

12-COMPARTMENT BIOPSY CONTAINER Half Prefilled with 10% Formalin

Made of polystyrene

This model will contain up to twelve biopsies. The screw on lid incorporates an O-ring in order to make it leakproof and protect its contents. The CoreDish™ conforms to OSHA directives. Compartments are numbered from 1 to 12.

Cat. #	For	Compartments	Prefilled	Qty/Pk	Qty/Cs
M970-D12	General purpose	12	YES	1	10
M971-D12	General purpose	12	NO	1	10

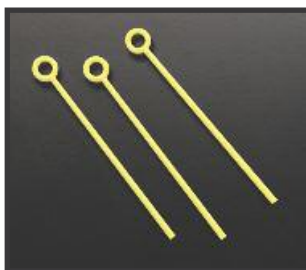

M975

The CorePicker™

Made of polystyrene

A practical tool and a great help to pick up and handle biopsies out of the CoreDish™. Packed in tamperproof resealable bags.

Cat. #	Length	Qty/Bag	Qty/Pk
M975	53 mm (2 1/8 in.)	25	125



M976

Shipping Box for CoreDish™

Made of cardboard

This sturdy and easy-to-assemble shipping box is most convenient for transporting or mailing the Simport® CoreDish.

Cat. #	Description	Qty/Pk
M976	Cardboard Box	10

CoreDish and CorePicker are trademarks of Simport® Scientific.

Simport®
Scientific

2588 Bernard-Pilon
Beloeil, Qc J3G 4S5 Canada

Phone: (450) 464-1723

Fax: (450) 464-3394

Email: info@simport.com

Website: www.simport.com

Distributed by: