



T360 OneHand™ Microtube Rack

Made of acetal polymer

A newly designed microtube rack that can be used all around the lab. Although one of the most attractive racks available today, it offers all the advantages required by the modern laboratory. Made of highly resistant acetal polymer, it will not shatter or stain in contact with most laboratory chemicals. No coating to worry about, which can chip, peel or rust in a water bath.

The OneHand™ Microtube Rack is compact, lightweight and stackable in order to save as much space as possible. This is why it is ideal for incubators, refrigerators, freezers, under lab hoods and on bench tops. Not only is it submersible but will also sink and maintain stability without tipping over.

The OneHand™ Microtube Rack is made of 2 tiers to facilitate the insertion and stability of microtubes. Now with one hand, you can easily unscrew just about any type of microcentrifuge tube with a screw-on closure. Thanks to an innovative locking system, the Simport self-standing Microwtubes® will securely lock in each well of the base tier and will not turn. All models of microtubes made by various manufacturers will lock in the upper tier, thanks to a series of teeth grasping the collar of the microtube. Convenient handles on each side of the rack will ensure a safe grip when carrying it around. Interlocking feet allow safe stacking. Available in 5 attractive colors. Individually wrapped.

Size: 293 mm x 115 mm x 39 mm H
(11 1/2 x 4 1/2 x 1 1/2 in. H)

ONE HAND OPENING AND CLOSING OF MICROTUBES




Tubes with ribs lock in place when engaged in serrated holes.



The locking base of the OneHand™ Rack locks self standing Simport Microwtubes in place.



Cat. #	Capacity	Color	Qty/Cs
T360-50B	50	Blue	10
T360-50G	50	Green	10
T360-50L	50	Lilac	10
T360-50O	50	Orange	10
T360-50Y	50	Yellow	10

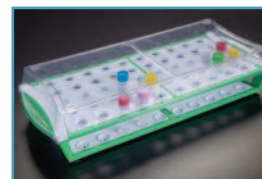

A transparent lid can be placed on the OneHand™, allowing a clear view of the contents.

Cat. #	Description	Color	Qty/Cs
T361-50	Cover	Clear	10



T350 MICREW TUBE® Storage Box

Made of polycarbonate

COLOR YOUR WORLD with this 100-place MICREW TUBE® Storage Box for tubes ranging from 0.5 ml. to 2 ml (except non self-standing 1.5 ml conical bottom tubes). Made of extra strong polycarbonate, this durable storage box is designed to be used at temperatures between -196 °C and +121 °C and is autoclavable.

A transparent cover allows the user to see the contents of the box, and is keyed to the base to prevent misalignment. Printed with a series of squares (numbered from 1 to 100), the surface accepts writing with markers, for better inventory control.

Samples can be classified more easily, thanks to a series of colored bases. A choice of four popular colors is available: blue, green, red and yellow.

Size: 133 mm x 133 mm x 52 mm H (5 1/4 x 5 1/4 x 2 1/16 in. H)

Cat. No.	For Micrewtubes	Color of base	Qty/Pk	Qty/Cs
T350-100B	0.5 to 2 ml	Blue	4	24
T350-100G	0.5 to 2 ml	Green	4	24
T350-100R	0.5 to 2 ml	Red	4	24
T350-100Y	0.5 to 2 ml	Yellow	4	24