



Simport[®]

Since 1975

Scientific inc.

Simport

Cassettes for the
Primera[®] cassette printers

In **QuickLoad[™] Stacks**



Available in QuickLoad Stacks or in Bulk.

www.simport.com


Specially designed for the automated and manual Primera Cassette printers.

Cassettes for Primera® Printers


Made of acetal

Disposable plastic cassettes hold tissue specimens during processing and embedding, as well as in storage. Molded from a special high-density polymer, these cassettes keep specimens safely submerged in liquid and are resistant to the chemical action of most histological solvents. **The efficient flow-through slots maximize fluid exchange and ensure proper drainage.** As required they can be opened or closed and they always relock securely without danger of specimen loss. Available in 11 colors.


HISTOSETTE II®

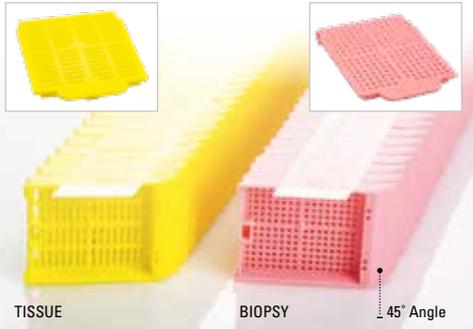

UNISETTE™ II

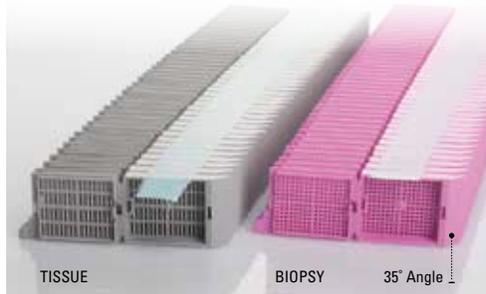

SLIMSETTE™ II

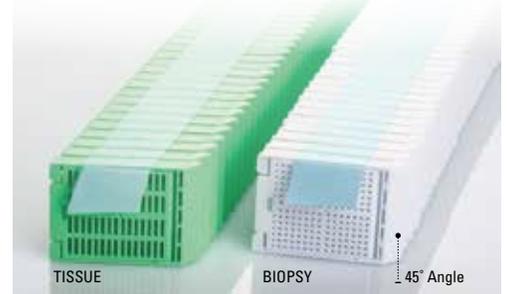

Automated feed printer		Manual feed printer			
50 stacks per case Case 2000		Individual Cassette 2 dispenser boxes of 500 cassettes and 1 dispenser box of 1000 covers / Case 1000			
TISSUE Product #	BIOPSY Product #	TISSUE Product #	BIOPSY Product #	Color	
M385-2T	M386-2T	M385-2	M386-2	<input type="checkbox"/>	White
M385-3T	M386-3T	M385-3	M386-3	<input type="checkbox"/>	Pink
M385-4T	M386-4T	M385-4	M386-4	<input type="checkbox"/>	Green
M385-5T	M386-5T	M385-5	M386-5	<input type="checkbox"/>	Yellow
M385-6T	M386-6T	M385-6	M386-6	<input type="checkbox"/>	Blue
M385-7T	M386-7T	M385-7	M386-7	<input type="checkbox"/>	Peach
M385-8T	M386-8T	M385-8	M386-8	<input type="checkbox"/>	Tan
M385-9T	M386-9T	M385-9	M386-9	<input type="checkbox"/>	Gray
M385-10T	M386-10T	M385-10	M386-10	<input type="checkbox"/>	Lilac
M385-11T	M386-11T	M385-11	M386-11	<input type="checkbox"/>	Orange
M385-12T	M386-12T	M385-12	M386-12	<input type="checkbox"/>	Aqua

Automated feed printer		Manual feed printer			
25 stacks per case Case 1000		Individual Cassette 3 dispenser boxes of 500 cassettes / Case 1500			
TISSUE Product #	BIOPSY Product #	TISSUE Product #	BIOPSY Product #	Color	
M405-2T	M406-2T	M405-2	M406-2	<input type="checkbox"/>	White
M405-3T	M406-3T	M405-3	M406-3	<input type="checkbox"/>	Pink
M405-4T	M406-4T	M405-4	M406-4	<input type="checkbox"/>	Green
M405-5T	M406-5T	M405-5	M406-5	<input type="checkbox"/>	Yellow
M405-6T	M406-6T	M405-6	M406-6	<input type="checkbox"/>	Blue
M405-7T	M406-7T	M405-7	M406-7	<input type="checkbox"/>	Peach
M405-8T	M406-8T	M405-8	M406-8	<input type="checkbox"/>	Tan
M405-9T	M406-9T	M405-9	M406-9	<input type="checkbox"/>	Gray
M405-10T	M406-10T	M405-10	M406-10	<input type="checkbox"/>	Lilac
M405-11T	M406-11T	M405-11	M406-11	<input type="checkbox"/>	Orange
M405-12T	M406-12T	M405-12	M406-12	<input type="checkbox"/>	Aqua

Automated feed printer		Manual feed printer			
50 stacks per case Case 2000		Individual Cassette 3 dispenser boxes of 500 cassettes / Case 1500			
TISSUE Product #	BIOPSY Product #	TISSUE Product #	BIOPSY Product #	Color	
M409-2T	M410-2T	M409-2	M410-2	<input type="checkbox"/>	White
M409-3T	M410-3T	M409-3	M410-3	<input type="checkbox"/>	Pink
M409-4T	M410-4T	M409-4	M410-4	<input type="checkbox"/>	Green
M409-5T	M410-5T	M409-5	M410-5	<input type="checkbox"/>	Yellow
M409-6T	M410-6T	M409-6	M410-6	<input type="checkbox"/>	Blue
M409-7T	M410-7T	M409-7	M410-7	<input type="checkbox"/>	Peach
M409-8T	M410-8T	M409-8	M410-8	<input type="checkbox"/>	Tan
M409-9T	M410-9T	M409-9	M410-9	<input type="checkbox"/>	Gray
M409-10T	M410-10T	M409-10	M410-10	<input type="checkbox"/>	Lilac
M409-11T	M410-11T	M409-11	M410-11	<input type="checkbox"/>	Orange
M409-12T	M410-12T	M409-12	M410-12	<input type="checkbox"/>	Aqua

M407 MICROMESH™ for Primera® Printers

Made of acetal


The MICROMESH™ offers 1676 square openings (0.38 mm) allowing for a greatly improved fluid exchange without having to use biopsy pads. Large anterior and posterior slots in both cassette and cover ensure that the MICROMESH™ Biopsy Cassette will sink rapidly. A large square compartment with a side measuring 27mm is perfect even for needle biopsies. The cover does not protrude above the cassette, a great space saving feature.

Molded from a special high density polymer, these cassettes keep specimens safely submerged in liquid and are totally resistant to the chemical action of most of histological solvents. The MICROMESH™ mesh ensures efficient fluid exchange and drainage. The one-piece integral lid eliminates the need for separate steel lids.

They can be opened and closed and always relock securely without danger of specimen loss. Anterior writing area is at a 45° angle.

Anatomy of Micromesh™ Cassettes

A recessed cover is a great space saving feature, allowing more cassettes to be stacked in automatic labeling machines and storage drawers

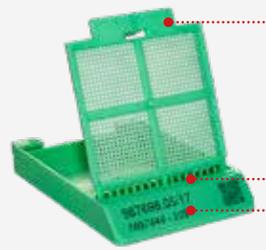

- Cover pre-installed on cassette
- Large anterior and posterior slots in both cassette and cover ensure that the Micromesh™ biopsy cassette will be submerged rapidly
- Anterior writing area is at a 45° angle submerged rapidly



Individual Cassette for manual feed printer

QuickLoad™ Stacks for Robotic feed printer



M407/M407T MICROMESH™		
Individual Cassette for manual feed printer 4 dispenser boxes of 250 / Case 1000	QuickLoad™ Stacks for robotic feed printer 50 stacks per case / Case 2000	Color
M407-2	M407-2T	<input type="checkbox"/> White
M407-3	M407-3T	<input type="checkbox"/> Pink
M407-4	M407-4T	<input type="checkbox"/> Green
M407-5	M407-5T	<input type="checkbox"/> Yellow
M407-6	M407-6T	<input type="checkbox"/> Blue
M407-7	M407-7T	<input type="checkbox"/> Peach
M407-8	M407-8T	<input type="checkbox"/> Tan
M407-9	M407-9T	<input type="checkbox"/> Gray
M407-10	M407-10T	<input type="checkbox"/> Lilac
M407-11	M407-11T	<input type="checkbox"/> Orange
M407-12	M407-12T	<input type="checkbox"/> Aqua

Primera® is a registered trademark of Primera Technology Incorporated



2588 Bernard-Pilon
Beloil, Quebec
J3G 4S5
Canada

T - 450 464-1723
F - 450 464-3394
info@simport.com
www.simport.com