


# Simport®

Canadian and Family Owned Since 1975

## Amplate™ Thin Wall PCR Plate


### NOW with Black Alphanumeric Grid!

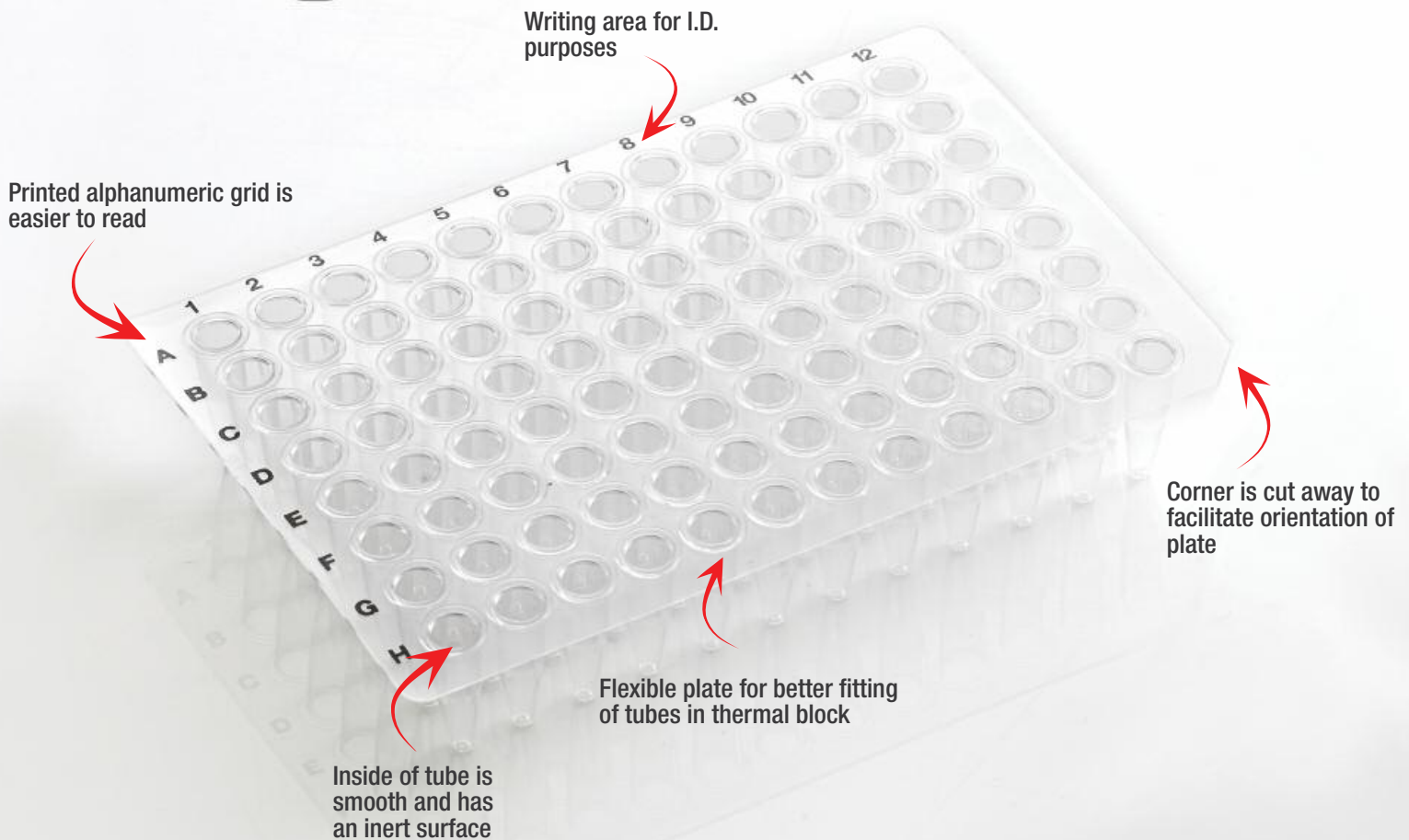

[www.simport.com](http://www.simport.com)

# Amplate™ Thin Wall PCR Plate

## NOW with Black Alphanumeric Grid!

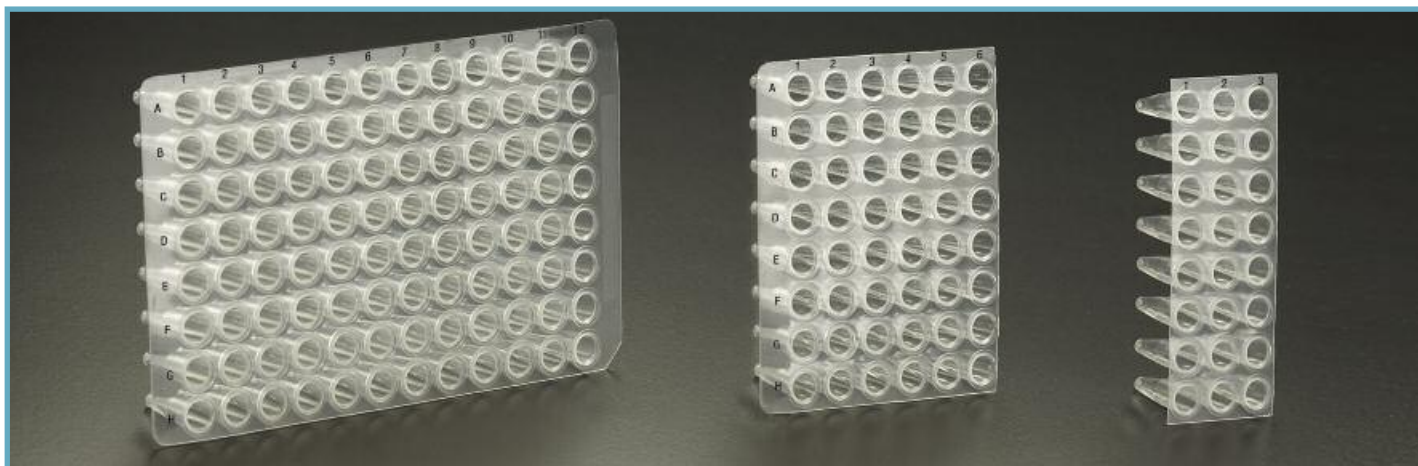

### T323 Amplate™ THIN WALL PCR PLATES

Made of polypropylene

Plates are certified RNase, DNase, Pyrogen and DNA-free

These 96-well, 48-well and 24-well PCR plates are thin-walled and designed for rapid thermal transfer. Each well has a capacity of 0.2 ml. They are precision-molded to ensure well-to-well and plate-to-plate uniformity. The inside of the wells is smooth and has an inert surface on which enzymes and nucleic acids do not bind.

Sealing can be performed with SecureSeal™ Thermal sealing film and foil (T329 Series), and Amplate™ Mat (T329-10). Suitable to be used with most 96-well thermal cyclers.

On these plates, a printed alphanumeric grid helps sample identification. To facilitate orientation, the bottom right corner of the plate is cut away. The AMPLATE™ is easy to seal since no cylindrical walls extend above the plate. More economical than using single tubes, it is available colorless and in four different colors. Packed in tamperproof resealable bags of 10 plates.

Cat. #	Color	Qty/Bag	Qty/Cs	Cat. #	Cat. #	Color	Qty/Bag	Qty/Cs
T323-96N	Natural	10	100	T323-48N	T323-24N	Natural	10	50
T323-96B	Blue	10	100	T323-48B	T323-24B	Blue	10	50
T323-96G	Green	10	100	T323-48G	T323-24G	Green	10	50
T323-96R	Red	10	100	T323-48R	T323-24R	Red	10	50
T323-96Y	Yellow	10	100	T323-48Y	T323-24Y	Yellow	10	50



### T323-96LP

#### Low Profile Amplate™ THIN WALL PCR PLATES

Made of polypropylene

Low profile 96-well PCR plates are also available. They are similar to the regular Simport® AMPLATE™ Series. However, each of the 96 wells has a smaller volume (only 100 µl) and thereby reduces the dead space between sample and cover.

Cat. #	Color	Qty/Bag	Qty/Cs
T323-96LPN	Natural	10	100
T323-96LPB*	Blue	10	100
T323-96LPG*	Green	10	100
T323-96LPR*	Red	10	100
T323-96LPY*	Yellow	10	100

\*Minimum quantity applicable. Please contact one of our customer service agents for further details.

Amplate is trademark of Simport® Scientific.

**Simport**<sup>®</sup>  
Scientific

Phone: (450) 464-1723

Fax: (450) 464-3394

Email: info@simport.com

Website: www.simport.com

2588 Bernard-Pilon  
Beloil, Qc J3G 4S5 Canada

Distributed by: